


Understanding Instech Sterile Pouches

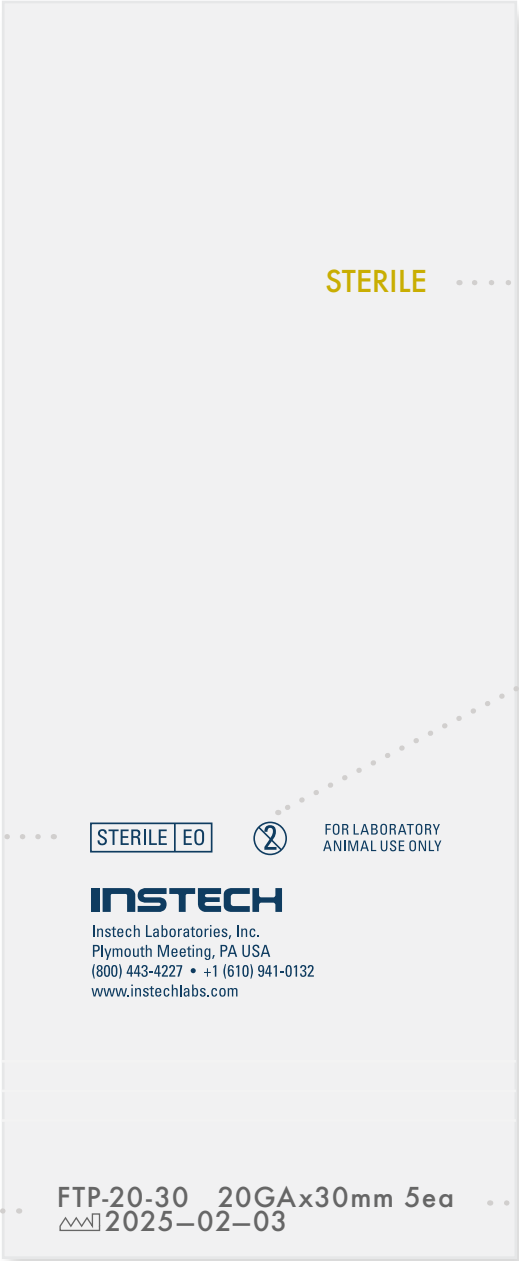

STERILE

Special ink turns from clear to brown when exposed to ethylene oxide gas.

ISO symbol indicating product has been sterilized by ethylene oxide gas.

STERILE EO



FOR LABORATORY ANIMAL USE ONLY

ISO symbol for "do not reuse." Many Instech products cannot be properly cleaned after use and most plastic parts will be damaged by autoclaving.

INSTECH

Instech Laboratories, Inc.
Plymouth Meeting, PA USA
(800) 443-4227 • +1 (610) 941-0132
www.instechlabs.com

FTP-20-30 20GAx30mm 5ea
2025-02-03

PART NUMBER: Description

Date of manufacture (YYYY-MM-DD). Instech's policy is to label its sterile products with the date of manufacture, not an expiration date. The actual time period after which sterility can become compromised depends on storage conditions, and the risk that non-sterility poses varies depending on the application. Therefore, Instech leaves the decision on when a product is no longer fit for use to its customers. A two-year shelf-life is typical for human medical devices but is fairly conservative; the manufacturers of the sterilizers and the pouches that Instech uses suggest a shelf-life of five years. This date of manufacture combined with the part number also serves as the lot number for our record keeping.